
Product Manual


Free Glycerol Assay Kit (Fluorometric)

Catalog Number

STA-399

100 assays

FOR RESEARCH USE ONLY
Not for use in diagnostic procedures

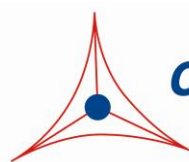

CELL BIOLABS, INC.

Creating Solutions for Life Science Research

Introduction

Glycerol is the backbone of Triglycerides (TAG). Triglycerides are a type of lipid in the blood, serving as an energy source and playing a key role in metabolism. Triglycerides are the digestive end product of breaking down dietary fats. Any extra carbohydrates and fats that are not immediately used are chemically converted into triglycerides. In the intestines, secreted enzyme lipases hydrolyse the triglyceride ester bond, yielding glycerol and free fatty acids in a process called lipolysis. Enterocytes then absorb and repackage the fragments with cholesterol to form chylomicrons, a major lipoprotein transport particle. In the liver, hepatic lipases also break down triglycerides to assemble another lipoprotein particle (VLDL) from triglycerides, cholesterol, and apolipoproteins.

Cell Biolabs' Free Glycerol Assay Kit measures free, endogenous glycerol by a coupled enzymatic reaction system. The glycerol is phosphorylated and oxidized, producing hydrogen peroxide which reacts with the kit's Fluorometric Probe (Ex. 530-560 nm/Em. 585-595 nm).

The Free Glycerol Assay Kit is a simple, fluorometric assay that quantitatively measures the amount of glycerol in plasma or serum in a 96-well microtiter plate format. Each kit provides sufficient reagents to perform up to 100 assays, including blanks, glycerol standards and unknown samples. The kit contains a glycerol standard and has a detection sensitivity limit of $\sim 1 \mu\text{M}$ (0.009 mg/dL).

Related Products

1. STA-369: OxiSelect™ Human Oxidized LDL ELISA Kit (MDA-LDL Quantitation)
2. STA-375: Uric Acid/UriCase Assay Kit
3. STA-378: Creatinine Assay Kit
4. STA-390: Total Cholesterol Assay Kit
5. STA-391: HDL and LDL/VLDL Cholesterol Assay Kit
6. STA-394: HDL Cholesterol Assay Kit
7. STA-398: Free Glycerol Assay Kit (Colorimetric)

Kit Components

1. Glycerol Standard (1 M) (Part No. 239801): One 200 μL vial of a 1 M glycerol solution.
2. 10X Assay Buffer (Part No. 239802): One 1.5 mL vial.
3. 5X Enzyme Mixture (Part No. 239803): Four 525 μL vials.
4. 200X Fluorometric Probe (Part No. 239901): One 55 μL amber vial.

Materials Not Supplied

1. Standard 96-well fluorescence black microtiter plate
2. 10 μL to 1000 μL adjustable single channel micropipettes with disposable tips
3. 50 μL to 300 μL adjustable multichannel micropipette with disposable tips
4. Multichannel micropipette reservoir

- Fluorescence microplate reader capable of reading excitation in the 530-560 nm range and emission in the 585-595 nm range.

Storage

Store entire kit at -80°C. Avoid multiple freeze/thaws by aliquoting. The Fluorometric Probe is light sensitive and should be maintained in amber tubes.

Preparation of Reagents

- Glycerol Standard, 10X Assay Buffer, and 5X Enzyme Mixture should be thawed/maintained at 4°C during assay preparation. All are stable for 1 week at 4°C. For longer term storage, each should be aliquoted and frozen at -80°C to avoid multiple freeze/thaws.
- 200X Fluorometric Probe should be thawed/maintained at room temperature during assay preparation. Any unused material should be aliquoted and frozen at -80°C to avoid multiple freeze/thaws.

Preparation of Glycerol Standard

- To prepare the glycerol standards, first perform a 1:100 dilution of the stock 1 M Glycerol Standard in deionized water. Prepare only enough for immediate use (e.g., Add 10 µL of 1 M Glycerol Standard to 990 µL deionized water). This solution has a concentration of 10 mM. Use this 10 mM glycerol solution to prepare standards in the concentration range of 0 µM – 100 µM by further diluting in deionized water (e.g., Add 10 µL of 10 mM glycerol solution to 990 µL deionized water - see Table 1 below). Glycerol diluted solutions and standards should be prepared fresh.

Standard Tubes	10 mM Glycerol Standard (µL)	Deionized Water (µL)	Final Glycerol Standard (µM)	Final Glycerol Standard (mg/dL)
1	10	990	100	0.92
2	250 of Tube #1	250	50	0.46
3	250 of Tube #2	250	25	0.23
4	250 of Tube #3	250	12.5	0.12
5	250 of Tube #4	250	6.25	0.06
6	250 of Tube #5	250	3.12	0.03
7	250 of Tube #6	250	1.56	0.015
8	0	250	0	0

Table 1. Preparation of Glycerol Standards

Preparation of Samples

- Plasma: Collect blood with an anticoagulant such as heparin, citrate or EDTA and mix by inversion. Centrifuge the blood at 1000 x g at 4°C for 10 minutes. Collect plasma supernatant without disturbing the white buffy layer. Sample should be tested immediately or frozen at -80°C for storage. Plasma may need to be diluted (up to 1:4 in PBS) before assaying. Normal human plasma has a glycerol concentration in the range of 0.12-0.61 mg/dL.

- Serum: Collect blood in a tube with no anticoagulant. Allow the blood to clot at room temperature for 30 minutes. Centrifuge at 2500 x g for 20 minutes. Remove the yellow serum supernatant without disturbing the white buffy layer. Samples should be tested immediately or frozen at -80°C for storage. Serum may need to be diluted (up to 1:4 in PBS) before assaying. Normal human serum has a glycerol concentration in the range of 0.4-1.2 mg/dL.

Assay Protocol

Each glycerol standard and sample should be assayed in duplicate or triplicate. A freshly prepared standard curve should be used each time the assay is performed.

1. Add 10 µL of the diluted glycerol standards or samples to the 96-well fluorescence microtiter plate.
2. **Maintain all components/mixtures at 4°C.** According to Table 2 (below), prepare the desired volume of Reaction Mixture (based on the # of tests) in the following sequence:
 - a. In a tube, add the appropriate volume of deionized water.
 - b. To the water add the corresponding volume of 10X Assay Buffer. Mix well.
 - c. Next, add the corresponding volume of 5X Enzyme Mixture.
 - d. Finally, add the corresponding volume of 200X Fluorometric Probe. Mix well and immediately use.

Note: Reaction Mixture will appear slightly pink in color. This is normal background and should be subtracted from all absorbance values.

Deionized Water (mL)	10X Assay Buffer (mL)	5X Enzyme Mixture (mL)	200X Fluorometric Probe (µL)	Total Volume of Reaction Mixture (mL)	# of Tests in 96-well Plate (90 µL/test)
5.950	1	2	50	9	100
2.975	0.5	1	25	4.5	50
1.190	0.2	0.4	10	1.8	20

Table 2. Preparation of Reaction Mixture

3. Transfer 90 µL of the above Reaction Mixture to each well (already containing 10 µL of glycerol standard or sample).
4. Cover the plate wells to protect the reaction from light.
5. Incubate at room temperature for 15 minutes on an orbital shaker.
6. Read the plate with a fluorescence microplate reader equipped for excitation in the 530-560 nm range and for emission in the 585-595 nm range.
7. Calculate the concentration of glycerol within samples by comparing the sample fluorescence to the standard curve. Negative controls (without glycerol) should be subtracted.

Example of Results

The following figures demonstrate typical Free Glycerol Assay Kit results. One should use the data below for reference only. This data should not be used to interpret actual results.

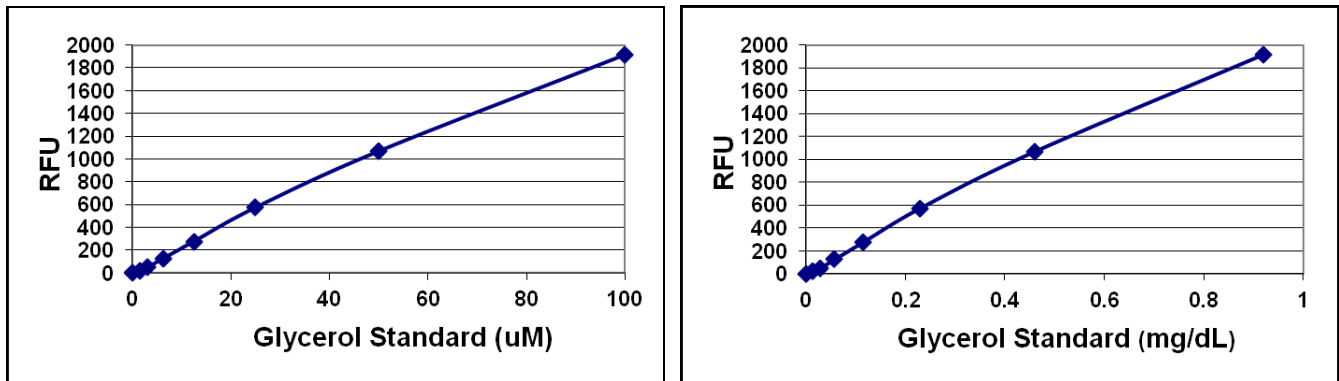


Figure 1: Free Glycerol Assay Standard Curve. Glycerol standard curve was performed according to the Assay Protocol, data plotted in μM (left) and mg/dL (right).

References

1. Bortz, W.M., Paul, P., Hff, A.C., and Holmes, W.H. (1972) *J. Clin. Invest.* **51**, 1537-1546.
2. Arner, P. (1995) *Internation J. of Obesity and Related Metabolic Disorders* **19(7)**, 435-442.
3. Ou, C-N and Frawley, V.L. (1985) *Clin. Biochem.* **18(1)**, 37-39.
4. Bernert, J.T., Bell, C.J., McGuffey, J.E., and Waymack, P.P. (1992) *J. Chromatogr.* **578**, 1-7.

Warranty

These products are warranted to perform as described in their labeling and in Cell Biolabs literature when used in accordance with their instructions. THERE ARE NO WARRANTIES THAT EXTEND BEYOND THIS EXPRESSED WARRANTY AND CELL BIOLABS DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR PARTICULAR PURPOSE. CELL BIOLABS's sole obligation and purchaser's exclusive remedy for breach of this warranty shall be, at the option of CELL BIOLABS, to repair or replace the products. In no event shall CELL BIOLABS be liable for any proximate, incidental or consequential damages in connection with the products.

Contact Information

Cell Biolabs, Inc.
 5628 Copley Drive
 San Diego, CA 92111
 Worldwide: +1 858 271-6500
 USA Toll-Free: 1-888-CBL-0505
 E-mail: tech@cellbiolabs.com
www.cellbiolabs.com

©2012-2023: Cell Biolabs, Inc. - All rights reserved. No part of these works may be reproduced in any form without permissions in writing.