# TDP-43/TARDBP Rabbit pAb

Catalog No.: A13405 1 Publications



## **Basic Information**

#### **Observed MW**

45kDa

#### **Calculated MW**

45kDa

#### Category

Polyclonal Antibody

#### **Applications**

WB,IHC-P,IF/ICC,ELISA

## **Cross-Reactivity**

Human, Mouse, Rat

## **Background**

HIV-1, the causative agent of acquired immunodeficiency syndrome (AIDS), contains an RNA genome that produces a chromosomally integrated DNA during the replicative cycle. Activation of HIV-1 gene expression by the transactivator Tat is dependent on an RNA regulatory element (TAR) located downstream of the transcription initiation site. The protein encoded by this gene is a transcriptional repressor that binds to chromosomally integrated TAR DNA and represses HIV-1 transcription. In addition, this protein regulates alternate splicing of the CFTR gene. A similar pseudogene is present on chromosome 20.

## **Recommended Dilutions**

**WB** 1:500 - 1:2000

IHC-P 1:50 - 1:200

**IF/ICC** 1:50 - 1:200

**ELISA** Recommended starting

concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

## Contact

www.abclonal.com

# **Immunogen Information**

**Gene ID**23435

Swiss Prot
Q13148

### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 1-100 of human TDP-43/TARDBP(NP\_031401.1).

### **Synonyms**

ALS10; TDP-43; TDP-43/TARDBP


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

## Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

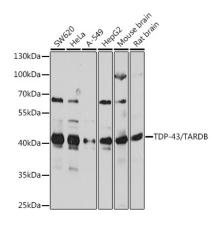

Western blot analysis of various lysates using TDP-43/TARDB Rabbit pAb (A13405) at 1:1000 dilution.

Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000


dilution.

Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 10s.

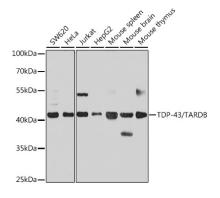

Western blot analysis of various lysates using TDP-43/TARDB Rabbit pAb (A13405) at 1:1000 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 30s.

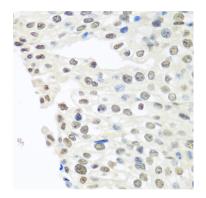

Immunohistochemistry analysis of paraffin-embedded Human prostate cancer using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

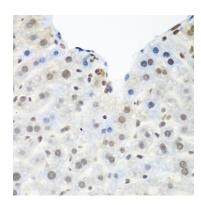


Immunohistochemistry analysis of paraffin-embedded Mouse liver using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

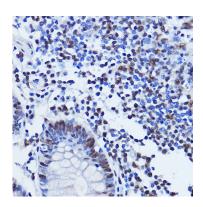


Immunohistochemistry analysis of paraffin-embedded Human colon using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:500 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.


## **Validation Data**

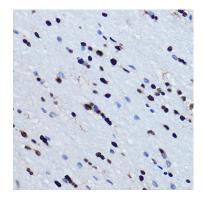


Immunohistochemistry analysis of paraffin-embedded Mouse brain using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:500 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

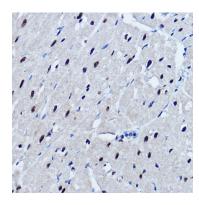


Immunohistochemistry analysis of paraffin-embedded Mouse heart using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:500 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

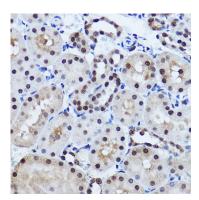


Immunohistochemistry analysis of paraffin-embedded Mouse kidney using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:500 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

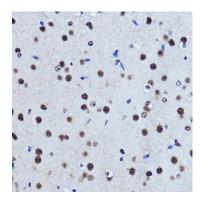


Immunohistochemistry analysis of paraffin-embedded Rat brain using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:500 (40x lens). Microwave antigen retrieval performed with 0.01M PBS Buffer (pH 7.2) prior to IHC staining.

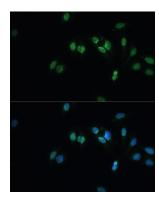


Immunofluorescence analysis of U-2 OS cells using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

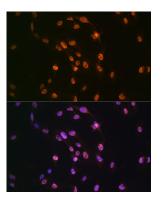

Immunofluorescence analysis of C6 cells using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

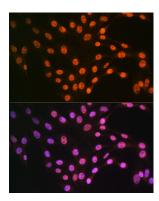

Immunofluorescence analysis of NIH-3T3 cells using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

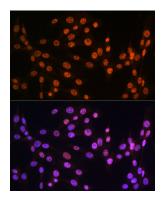

Immunofluorescence analysis of NIH-3T3 cells using TDP-43/TARDB Rabbit pAb (A13405) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.