

Neuroigin 1 Antibody
Neuroigin 1 Antibody, Clone S97A-31
Catalog # ASM10295

Specification

Neuroigin 1 Antibody - Product Information

| | |
|-------------------|-----------------------------|
| Application | WB, IHC, ICC |
| Primary Accession | Q62765 |
| Other Accession | NP_446320.1 |
| Host | Mouse |
| Isotype | IgG1 |
| Reactivity | Human, Mouse, Rat |
| Clonality | Monoclonal |

Description

Mouse Anti-Rat Neuroigin 1 Monoclonal IgG1

Target/Specificity

Detects ~120kDa. Does not cross-react with other Neuroiginins.

Other Names

NLG 1 Antibody, KIAA1070 Antibody, MGC45115 Antibody, Neuroigin-1 Antibody, NL1 Antibody, NLG1 Antibody, Nlgn1 Antibody

Immunogen

Fusion protein amino acids 718-843 (cytoplasmic C-terminus) of rat Neuroigin-1. Mouse: 99% identity (125/126 amino acids identical). Human: 99% identity (125/126 amino acids identical) >40% identity with Neuroigin-2 and -3.

Purification

Protein G Purified

Storage **-20°C**

Storage Buffer

PBS pH 7.4, 50% glycerol, 0.1% sodium azide

Shipping Temperature **Blue Ice or 4°C**

Certificate of Analysis

2 µg/ml of SMC-463 was sufficient for detection of Neuroigin-1 in 20 µg of rat brain lysate by colorimetric immunoblot analysis using Goat anti-mouse IgG:HRP as the secondary antibody.

Cellular Localization

Cell Membrane | Cell Junction | Synapse | Postsynaptic Cell Membrane | Postsynaptic Density


Neuroigin 1 Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.

- [Western Blot](#)
- [Blocking Peptides](#)

- [Dot Blot](#)
- [Immunohistochemistry](#)
- [Immunofluorescence](#)
- [Immunoprecipitation](#)
- [Flow Cytometry](#)
- [Cell Culture](#)

Neuroigin 1 Antibody - Images

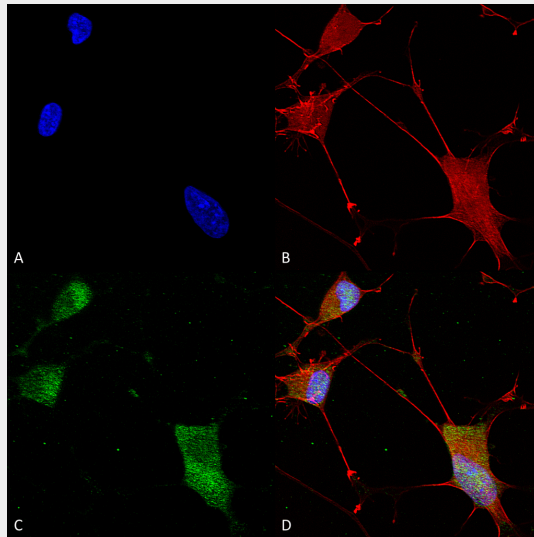

Immunocytochemistry/Immunofluorescence analysis using Mouse Anti-Neuroigin 1 Monoclonal Antibody, Clone N97A/31 (ASM10295). Tissue: Neuroblastoma cells (SH-SY5Y). Species: Human. Fixation: 4% PFA for 15 min. Primary Antibody: Mouse Anti-Neuroigin 1 Monoclonal Antibody (ASM10295) at 1:100 for overnight at 4°C with slow rocking. Secondary Antibody: AlexaFluor 488 at 1:1000 for 1 hour at RT. Counterstain: Phalloidin-iFluor 647 (red) F-Actin stain; Hoechst (blue) nuclear stain at 1:800, 1.6mM for 20 min at RT. (A) Hoechst (blue) nuclear stain. (B) Phalloidin-iFluor 647 (red) F-Actin stain. (C) Neuroigin 1 Antibody (D) Composite.

Neuroigin 1 Antibody - Background

Neuroigin-1 is a neuronal cell surface protein belonging to the type-B carboxylesterase/lipase family. It is a necessary component in the maturation of excitatory synapses for their normal, functional development, and is necessary to the regulation of synaptic plasticity and the development of long-term memory within the adult amygdala in mammals. It is believed to participate in cell-cell-interaction through the assembly of intracellular junction by the binding of beta-neurexins, and may also be a factor in the maintenance and assembly of synaptic junctions. It is also thought to have involvement in excitatory synaptic specification. Within brain tissue, Neuroigin-1 is primarily observed in neurons and spinal cord.