

A1443

Leader in Biomolecular Solutions for Life Science



CYP11B2 Rabbit pAb

Catalog No.: A1443

Basic Information

Observed MW

50kDa

Calculated MW

58kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,IF/ICC,ELISA

Cross-Reactivity

Human,Mouse,Rat

Background

This gene encodes a member of the cytochrome P450 superfamily of enzymes. The cytochrome P450 proteins are monooxygenases which catalyze many reactions involved in drug metabolism and synthesis of cholesterol, steroids and other lipids. This protein localizes to the mitochondrial inner membrane. The enzyme has steroid 18-hydroxylase activity to synthesize aldosterone and 18-oxocortisol as well as steroid 11 beta-hydroxylase activity. Mutations in this gene cause corticosterone methyl oxidase deficiency.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:100

IF/ICC 1:50 - 1:100

Immunogen Information

Gene ID

1585

Swiss Prot

P19099

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 234-503 of human CYP11B2 (NP_000489.3).

Synonyms

CPN2; ALDOS; CYP11B; CYP11BL; CYPXIB2; P450C18; P-450C18; P450aldo; CYP11B2

Contact

 www.abclonal.com

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

Validation Data

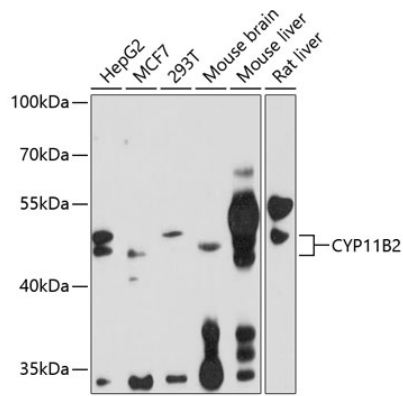

Western blot analysis of various lysates using CYP11B2 Rabbit pAb (A1443) at 1:1000 dilution.


Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Enhanced Kit (RM00021).


Exposure time: 20s.

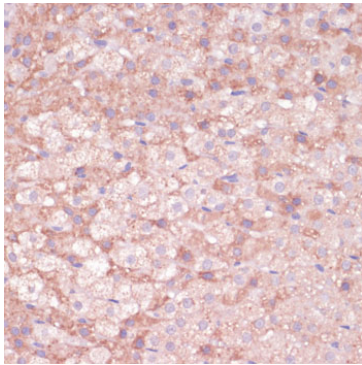

Immunohistochemistry analysis of paraffin-embedded rat adrenal gland using CYP11B2 Rabbit pAb (A1443) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

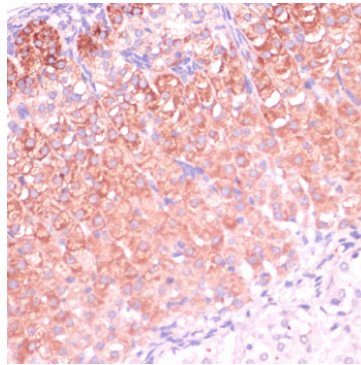

Immunohistochemistry analysis of paraffin-embedded mouse adrenal gland using CYP11B2 Rabbit pAb (A1443) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

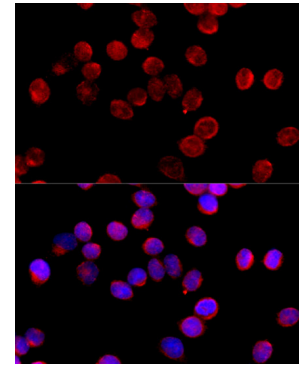

Immunofluorescence analysis of HeLa cells using CYP11B2 Rabbit pAb (A1443) at dilution of 1:100 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.