

# CD27 antibody

#### **Product Information**

Catalog No.:	FNab01453
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	$\geq$ 95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

#### Background

Receptor for CD70/CD27L. May play a role in survival of activated T-cells. May play a role in apoptosis through association with SIVA1.


### Immunogen information

Immunogen:	CD27 molecule
Synonyms:	CD27 antigen CD27L receptor T-cell activation antigen CD27 T14 Tumor necrosis factor receptor superfamily member 7 CD27 TNFRSF7
Observed MW:	50 kDa
Uniprot ID :	P26842

## Application

Reactivity:	Human			
Tested Application:	ELISA, WB, FC			
Recommended dilution: WB: 1:500-1:1000				
Image:				



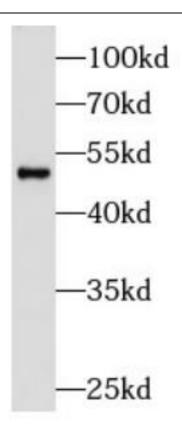

Jurkat cells were subjected to SDS PAGE followed by western blot with FNab01453(CD27 antibody) at dilution of 1:600