ABclonal www.abclonal.com

Histone H3 Rabbit pAb

Catalog No.: A2348 167 Publications

Basic Information

Observed MW

17kDa

Calculated MW

16kDa

Category

Polyclonal Antibody

Applications

WB,IHC-P,IF/ICC,IP,ChIP,ELISA

Cross-Reactivity

Human, Mouse, Rat, Other (Wide Range Predicted)

Background

Histones are basic nuclear proteins that are responsible for the nucleosome structure of the chromosomal fiber in eukaryotes. Nucleosomes consist of approximately 146 bp of DNA wrapped around a histone octamer composed of pairs of each of the four core histones (H2A, H2B, H3, and H4). The chromatin fiber is further compacted through the interaction of a linker histone, H1, with the DNA between the nucleosomes to form higher order chromatin structures. This gene is intronless and encodes a replication-dependent histone that is a member of the histone H3 family. Transcripts from this gene lack polyA tails; instead, they contain a palindromic termination element. This gene is located separately from the other H3 genes that are in the histone gene cluster on chromosome 6p22-p21.3.

Recommended Dilutions

WB 1:2000 - 1:10000

IHC-P 1:50 - 1:200

IF/ICC 1:50 - 1:200

IP 0.5μg-4μg antibody for 200μg-400μg extracts

of whole cells

ELISA Recommended starting

concentration is 1 µg/mL. Please optimize the concentration based on your specific

assay requirements.

ChIP 3µg antibody for

5μg-10μg of Chromatin

Contact

www.abclonal.com

Immunogen Information

 Gene ID
 Swiss Prot

 8290/8350
 Q16695/P68431

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 36-136 of human HIST3H3 (NP_003520.1).

Synonyms

H3t; H3.4; H3/g; H3FT; H3C16; HIST3H3; Histone H3


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

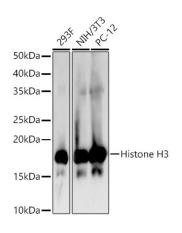

Western blot analysis of various lysates using Histone H3 Rabbit pAb (A2348) at 1:10000 dilution.


Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.

Lysates / proteins: 25 µg per lane.


Blocking buffer: 3 % nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).


Exposure time: 10s.

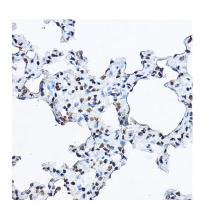


Immunohistochemistry analysis of paraffin-embedded Rat lung using Histone H3 Rabbit pAb (A2348) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

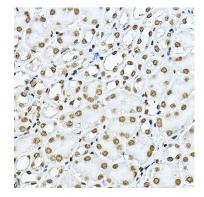


Immunohistochemistry analysis of paraffin-embedded Mouse kidney using Histone H3 Rabbit pAb (A2348) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

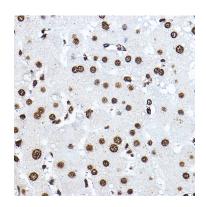


Immunohistochemistry analysis of paraffin-embedded Human liver using Histone H3 Rabbit pAb (A2348) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

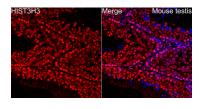

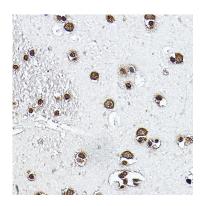

Immunohistochemistry analysis of paraffin-embedded Human brain using Histone H3 Rabbit pAb (A2348) at dilution of 1:100 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.

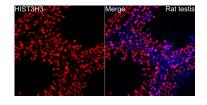


retrieval with 10 mM citrate buffer pH 6.0 before commencing with IF staining protocol.Immunofluorescence analysis of paraffin-embedded mouse testis using Histone H3 Rabbit pAb (A2348) at dilution of 1:200 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue:

Perform high pressure antigen


Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IF staining protocol.Immunofluorescence analysis of paraffin-embedded rat testis using Histone H3 Rabbit pAb (A2348) at dilution of 1:200 (40x

lens). Secondary antibody: Cy3conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.


DAPI for nuclear staining.

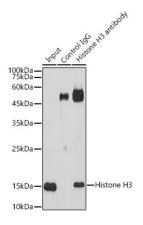

Immunoprecipitation analysis of 600 μ g extracts of HeLa cells using 10ug Histone H3 antibody (A2348). Western blot was performed from the immunoprecipitate using Histone H3 antibody (A2348) at a dilution of 1:30000.

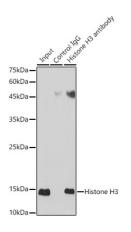

Immunoprecipitation analysis of 600 μ g extracts of HeLa cells using 10ug Histone H3 antibody (A2348). Western blot was performed from the immunoprecipitate using Histone H3 antibody (A2348) at a dilution of 1:30000.

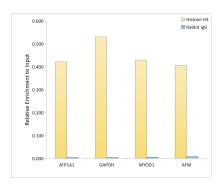

Chromatin immunoprecipitation analysis of extracts of 293T cells, using Histone H3 antibody (A2348) and rabbit IgG.The amount of immunoprecipitated DNA was checked by quantitative PCR. Histogram was constructed by the ratios of the immunoprecipitated DNA to the input.