

CFHR3 antibody

Product Information

Catalog No.:	FNab10775
Size:	100µg
Form:	liquid
Purification:	Immunogen affinity purified
Purity:	≥95% as determined by SDS-PAGE
Host:	Rabbit
Clonality:	polyclonal
Clone ID:	None
IsoType:	IgG
Storage:	PBS with 0.02% sodium azide and 50% glycerol pH 7.3, -20°C for 12 months(Avoid repeated freeze / thaw cycles.)

Background


The protein encoded by this gene is a secreted protein, which belongs to the complement factor H-related protein family. It binds to heparin, and may be involved in complement regulation. Mutations in this gene are associated with decreased risk of age-related macular degeneration, and with an increased risk of atypical hemolytic-uremic syndrome. Alternatively spliced transcript variants encoding different isoforms have been found for this gene.

Immunogen information


Immunogen:	CFHR3
Synonyms:	Complement factor H-related protein 3 (FHR-3) DOWN16 H factor-like protein 3 CFHR3 CFHL3 FHR3
Observed MW:	45 kDa
Uniprot ID :	Q02985

Application

Reactivity:	Human, Mouse
Tested Application:	ELISA, WB, IHC, IF
Recommended dilution:	WB: 1:500-1:2000; IHC: 1:50-1:500; IF: 1:50-1:100
Image:	

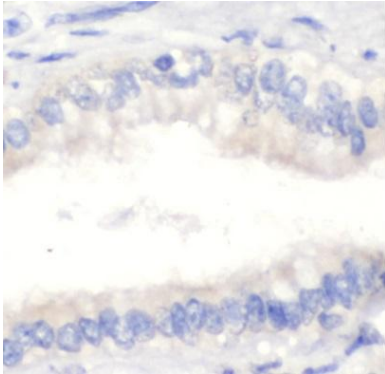

Immunohistochemistry of paraffin-embedded human kidney tissue slide using FNab10775(CFHR3 Antibody) at dilution of 1:100.

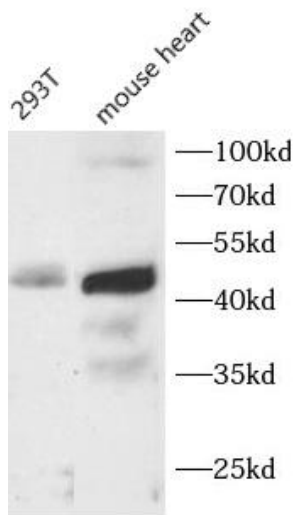

Various lysates were subjected to SDS PAGE followed by western blot with FNab10775(CFHR3 antibody) at dilution of 1:1000