

# Recombinant SARS-CoV-2 Spike S1+S2 ECD(S-ECD) Protein

Catalog No.: RP01283LQ Recombinant 3 Publications

# **Sequence Information**

**Species Gene ID Swiss Prot** HEK293 cells 43740568 P0DTC2

# Tags

C-His

### **Synonyms**

Envelope;SARS-CoV-2 Spike RBD (N501Y);Spike;Spike ECD;Spike RBD;Spike S1;Spike S2;Spike S2 ECD;S1-RBD protein;NCP-CoV RBD Protein;novel coronavirus RBD Protein;2019-nCoV RBD Protein;S glycoprotein Subunit1 RBD Protein

## **Product Information**

**Source** Purification
HEK293 cells > 95% by SDSPAGE.

### **Endotoxin**

< 0.1 EU/ $\mu g$  of the protein by LAL method.

# **Formulation**

Supplied as a 0.22 µm filtered solution in PBS, pH 7.4.Contact us for customized product form or formulation.

### Reconstitution

**Background** 

# **Basic Information**

### Description

Recombinant SARS-CoV-2 S1+S2 ECD(S-ECD) Protein is produced by HEK293 cells expression system. The target protein is expressed with sequence (Val11-Gln1208) of SARS-COV-2 S1+S2 ECD(S-ECD) (Accession #YP\_009724390.1) fused with a 6×His tag at the C-terminus.It is the wildtype sequence in pre-fusion state without any mutation.

# **Bio-Activity**

Measured by its binding ability in a functional ELISA. Immobilized SARS-CoV-2 S1+S2 ECD(S-ECD) at 2  $\mu$ g/mL (100  $\mu$ L/well) can bind recombinant Human ACE2 with a linear range of 0.15-3.72 ng/mL.

# Storage


Store at -70°C. This product is stable at  $\leq$  -70°C for up to 1 year from the date of receipt. For optimal storage, aliquot into smaller quantities after centrifugation and store at recommended temperature. Avoid repeated freeze-thaw cycles. Avoid repeated freeze/thaw cycles.

## **Contact**




www.abclonal.com

# **Validation Data**

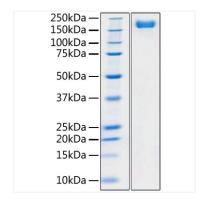

Recombinant SARS-COV-2 S1+S2 ECD(S-ECD) Protein with His tag was determined by SDS-PAGE with Coomassie Blue, showing a band at 180 kDa.

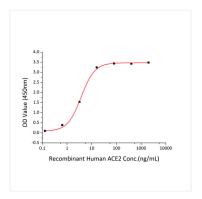

Immobilized SARS-CoV-2 S1+S2 ECD(S-ECD) at  $2\mu g/mL$  (100  $\mu L/well$ ) can bind recombinant Human ACE2 with a linear range of 0.15-3.72 ng/mL.